Foot Warmer with Dual Speed Massage

VM949J

Usage and Maintenance Instructions

This file is available to view and download as a PDF at www.aidapt.co.uk. Sight impaired customers can use a free PDF Reader (such as adobe.com/reader) to zoom in and increase the text size for improved readability.

Read all instructions before using this item:

- This product is supplied with a 12 Volt Adaptor. Make sure that the product is connected to an outlet having the same configuration as the plug.
- Always unplug the massager from the power source after use, and before cleaning.
- Do not reach for massager if it has fallen into water. Unplug immediately.
- An appliance should never be left unattended when plugged in.
- Unplug from outlet when not in use, before putting on or taking off parts or attachments.
- This product is intended for household use only. Do not use outdoors.
- Do not carry this appliance by supplied cord as a handle.
- Keep the cord away from heated surfaces.
- Close supervision is necessary when this appliance is used by, on or near children or those with limited dexterity or mobility.
- This unit is not a toy. Children should not use it or play with it.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use near eyes or other highly sensitive areas.
- Do not use heat for more than 30 minutes; do not exceed 30 minutes of use.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not provided with the unit.

BEFORE USE

- Carefully remove all packaging. Avoid using any knives or other sharp instruments as this may damage the product.
- Ensure that the product is dry and not folded or wrinkled before use.
- Check the product for any visible damage. If you see any damage or suspect a fault, do not use your product, but contact your supplier for support.

HEALTH CAUTION

- If you have any concerns regarding your health, consult your doctor before using this product.
- Individuals with pacemakers should not use this product.
- If you experience pain in a muscle or joint for a prolonged period of time. Discontinue use and consult your doctor. Persistent pain could be a symptom of a more serious condition.


INTENDED USE & FEATURES

The Foot Warmer with Dual Speed Massage has been designed to soothe cold, tired or aching feet and should be used for no more than 30 minutes at a time.

Features include:

- Soothing heat for added relief
- Two powerful massage motors provide soothing relief for your feet
- Massage and heat functions can be selected together or separately
- High and low massage levels
- Soft, machine-washable fleece lining

REMOTE CONTROL FUNCTIONS

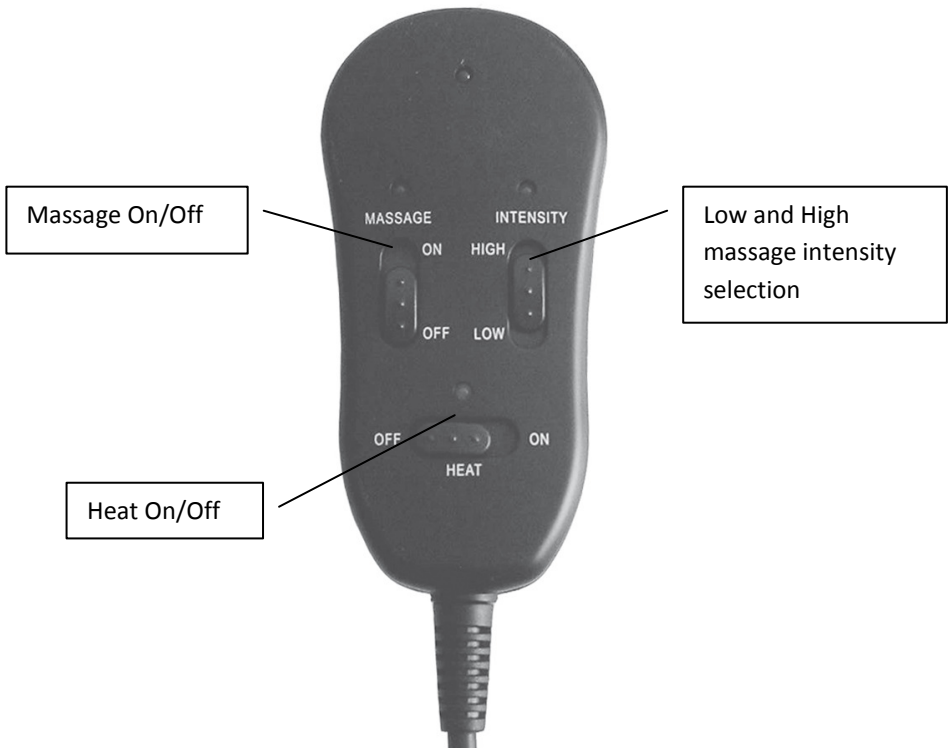CARE & MAINTENANCE

TO STORE

Place unit in its box or safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric. Do not wrap power cord around the unit. Do not hang the unit by power cord.

TO CLEAN THE LINING

1. Remove the fleece lining from the foot massager. The inner lining is washable, but do not attempt to machine wash outer lining.
2. Machine wash with a mild detergent. Use non-chlorine bleach only when necessary.
3. Do not dry clean.

TO CLEAN THE OUTER

1. Make sure the massager is unplugged.
2. Allow the unit to cool and the wipe down with a damp clean cloth.
3. Do not use any strong chemical abrasive cleaners.

IMPORTANT INFORMATION

The information given in this instruction booklet must not be taken as forming part of or establishing any contractual or other commitment by Aidapt Bathrooms Limited, Aidapt (Wales) Ltd or its agents or its subsidiaries and no warranty or representation concerning the information is given.

Please exercise common sense and do not take any unnecessary risks when using this product; as the user you must accept liability for safety when using the product.

Please do not hesitate to contact the person who issued this product to you or the manufacturer (detailed below) if you have any queries concerning the assembly/use of your product.

Aidapt Bathrooms Ltd, Lancots Lane, Sutton Oak, St Helens, WA9 3EX
Telephone: +44 (0) 1744 745 020 • Fax: +44 (0) 1744 745 001 • Web: www.aidapt.co.uk
Email: accounts@aidapt.co.uk • adaptations@aidapt.co.uk • sales@aidapt.co.uk

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Altai Germany (Altai Europe Ltd.)